

Oxygen Sensor OOM106

Use the advantages:

- Compliant with European MDD (CE certification)
- Meets ISO 80601-2-55
- Designed and manufactured according to EN ISO 13485
- Accurate and reliable response
- Resistant to N₂O
- Excellent signal stability
- High product quality
- Short delivery times
- Technical support
- · Made in Germany
- FDA cleared



From standard sensors to customized sensors

Experienced EnviteC engineers analyze customer requirements. This input is used for different standard and OEM applications, and ongoing support is provided right up to the final integration in the solution. EnviteC designs customized sensors characterized by a maximum possible degree of precision, for example with different signal levels or temperature compensation elements.

Intendend use

The EnviteC Medical Oxygen Sensors are intended as oxygensensing component of an oxygen analyzer that measures oxygen concentration in breathing gas mixtures in the following applications:


Sensing device for oxygen in

- control device of oxygen concentrators
- medical ventilators
- anaesthesia equipment
- incubators.

The use is limited to system monitoring. The sensors are not suited for breath by breath analysis of breath gases. Please refer to the Instructions for Use! If the sensor is intended to replace the original oxygen-sensing component of an oxygen analyzer, consult the EnviteC XRL Cross Reference List for selecting the appropriate sensor.



Mechanical drawing (All dimension in mm)

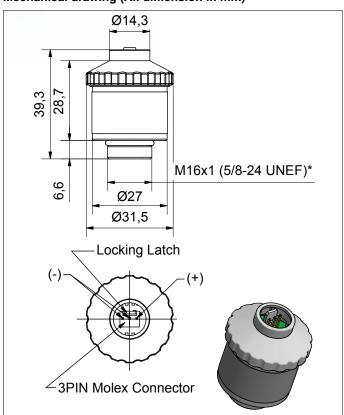

General tolerances ISO 2768-c

*Intermediate thread: Metric / Unified Extra Fine

Additional information

The Instructions for Use as well as the EnviteC XRL Cross Reference List are available under www.EnviteC.com and in the Apple App Store under EnviteC XRL as free download.

For more information please contact us! We look forward to assisting you either on the phone or in a personal talk.

Technical Specifications OOM106

Measurement range	0 % 100 % oxygen (at atmospheric pressure)
Nominal sensor lifetime	≥ 1 000 000 % volume oxygen hours
Output in ambient air	9 mV 13 mV
Electrical interface	3 pin (Molex® 22-11-1031)
Accuracy	meets ISO 80601-2-55 requirements
Repeatability	$<$ 1 % volume ${\rm O_2}$ at constant temperature and pressure
Linearity error	< 3 % relative
Response time	< 12 s to 90 % of final value
Zero offset voltage	< 200 μV in 100 % nitrogen, applied for 5 min
Cross interference	meets ISO 80601-2-55 requirements
Influence of humidity	-0.03 % rel. per % RH at 25 °C
Pressure range	0.6 bar 2 bar (ppO ₂ 0 1250 mbar O ₂)
Influence of pressure	proportional to change in oxygen partial pressure
Influence of mechanical shock	< 1 % relative after a fall from 1 m
Operating temperature	0 °C +50 °C
Temperature compensation	built-in NTC compensation
Effect of temperature compensation (steady	between +25 °C and +40 °C: 3 % relative error
state)	between 0 °C and +50 °C: 8 % relative error
Operating humidity	0 % 99 % RH non-condensing
Long term output drift	< 1 % volume oxygen per month
	typically < -15 % relative over lifetime
Storage temperature	-20 °C +50 °C
Recommended storage	+5 °C +15 °C
Recommended load	≥ 10 kOhms
Warm-up time	< 30 minutes, after replacement of sensor
Weight	approximately 28 grams
Part number	01-00-0091





All specifications are applicable at standard conditions: 1013 hPa, 25 $^{\circ}\text{C}$ dry ambient air

For suitable accessories and sensors please refer to the EnviteC Cross Reference List under www.EnviteC.com and in the Apple App Store unter EnviteC XRL as free download.

EnviteC-Wismar GmbH a Honeywell Company

Alter Holzhafen 18, 23966 Wismar, Germany

Phone: +49 (0)3841-360-1
Phone: +49 (0)3841-360-200
Fax: +49 (0)3841-360-222
Internet: www.envitec.com
Email: info@envitec.com

Doc. No. 001-33-Datasheet_OOM106-0 March 2016 Technical information is subject to change without notice! © 2016 Honeywell International Inc.

EnviteC by Honeywell reserves the right to make changes in product specifications and adjust its production at any time and without notice.

